

PREMIER MAGNESIA - GILES CHEMICAL

COMPANY PROCEDURE

USP Analytical Balance Uncertainty Measurement

Page : 1 of 2

Revision: 00 Effective Date: 02/21/2012

Author: Stephen Ballew Procedure Number: QA-LAB-44



Overview:

USP Section 41 states that "when substances are to be 'accurately weighed' for Assay the weighing is to be performed with a weighing device whose measurement uncertainty (random plus systematic error) does not exceed 0.1% of the reading." The formula for arriving at this measurement uncertainty value is three times the standard deviation, of at least ten replicate measurements, divided by the amount weighed. The resulting value shall not exceed 0.001.

Frequency:

An uncertainty measurement will be performed on the analytical balance by a quality analyst on a monthly basis.

Procedure:

- 1) Clean balance of any dust or powder using the brush supplied with the check weights, which are located in the marked box in the drawer below the analytical balance. After cleaning close all balance doors.
- 2) Press zero key and check that digital readout reads 0.00000 g within 4 seconds. Repeat once if digital readout fails to display 0.00000 g within 4 seconds.
- 3) Press Adjust.int function key. The balance will run through a calibration, zeroing check and temperature adjusting procedure. When it is finished it will signal 'Adjustment Done' on the screen along with the temperature. If 'Adjustment Done' does not appear within one minute, repeat the procedure once. If 'Adjustment Done' still does not appear report the problem to the laboratory supervisor or the person responsible for the balance and contact METTLER TOLEDO's service organization for advice.
- 4) Measure the uncertainty of the balance by carefully weighing the 5 g check weight 10 times. Always use the forceps provided to handle the calibration weights. Do NOT handle weights with fingers. Record the results.

Calculating Uncertainty:

(3 * Std. Dev. of ten measurements) ≤ 0.001 Amount Weighed

The data can then be entered into an excel spreadsheet. This can be used to calculate standard deviation and uncertainty. Record these results.

Evaluation:

If the uncertainty is above the maximum amount of 0.001 then repeat the test. If the uncertainty maximum is exceeded again, level the balance, perform adjustment with built-in or external adjustment weights and repeat the test. If the control limit is still exceeded, report the problem to the laboratory supervisor or the person responsible for the balance, mark the balance as "out of uncertainty limits" and contact METTLER TOLEDO's service organization for advice.



PREMIER MAGNESIA - GILES CHEMICAL

COMPANY PROCEDURE


USP Analytical Balance Uncertainty Measurement

Stephen

Ballew

Page : 2 of 2 Revi

Revision : 00 Effective Date : 02/21/2012

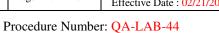



Revision Number	Revision Date	Effective Date	Revision Author	Quality Approval	Production Approval	Revision Description
00	02/21/2012	02/21/2012	Stephen Ballew	Deborah Durbin	Jason Bumgarner	New Document