

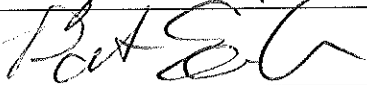

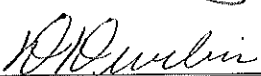
	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 1 of 10	

Approvals

Signing below indicates agreement that the protocol is ready for execution of the Installation and Operational Qualification for the Quincy Labs 20GC Drying Oven located at Giles Chemical QA Laboratory, 102 Commerce Street in Waynesville, NC.

Project Team Member	Functional Area	Signature	Date
Bryan Elchert	QA Laboratory		7/23/13
Patrick Owen	Engineering		7/24/13
Matt Haynes	Operations		7/24/13
Deborah Durbin	Quality		7/24/13

A final summary report that consists of results and conclusions based on the data collected after protocol execution will be written and approved. The executed protocol will be attached behind the report.

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.

**GILES CHEMICAL ~ PREMIER MAGNESIA****Validation Protocol**

Title: Drying Oven IQ/OQ Protocol

Number: L13-VAL-100-054

Owner: Bryan Elchert

Revision: 0

Effective Date: July 17, 2013



Page: 2 of 10



TABLE OF CONTENTS		Page #
APPROVAL PAGE		1
TABLE OF CONTENTS		2
I.	PURPOSE	3
II.	BACKGROUND	3
III.	OVERVIEW	3
IV.	SYSTEM DESCRIPTION	3
V.	SCOPE	3
VI.	ROLES AND RESPONSIBILITIES	3
VII.	TEST PROGRAM	4-5
A	INSTALLATION QUALIFICATION	4
B	OPERATIONAL QUALIFICATION	5
VIII.	CALIBRATION	5
ATTACHMENT I:	INSTALLATION QUALIFICATION	6
ATTACHMENT II:	OPERATIONAL QUALIFICATION	7
ATTACHMENT III:	PROTOCOL DEVIATION REPORT LOG	8
ATTACHMENT IV:	PROTOCOL DEVIATION REPORT	9
ATTACHMENT V:	SIGNATURE IDENTIFICATION LOG SHEET	10

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.

	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 3 of 10	

I. PURPOSE:

The purpose of this protocol is to certify with documented evidence that the Drying Oven is installed and functions as intended. This protocol sets forth the objectives, methodology, documentation, and test activities needed to complete the Installation Qualification (IQ) and Operational Qualification (OQ) for the Quincy Labs 20GC Drying Oven located at Giles Chemical QA Laboratory, 102 Commerce Street, Waynesville, NC.

II. BACKGROUND:

Giles Chemical, a division of Premier Magnesia (Giles), is dedicated to offering high quality Magnesium Sulfate products. To help achieve this goal, product samples sometimes must be efficiently dried for laboratory analysis for quality purposes.

III. OVERVIEW

No other departments or systems will be affected by the installation or use of this equipment.

The following tests will be performed in this qualification:

Installation Documentation – the serial number or asset tag number of the Drying Oven will be documented.

Utility Verification – the voltage to the Drying Oven will be documented and verified to be correct.

Control / Operation Verification – the controls will be verified.

IV. SYSTEM DESCRIPTION:

A. The system consists of a Quincy Labs 20GC Drying Oven.

B. Description of Operation

01. The Drying Oven is started by pressing the power button.

02. The Drying Oven heats up to desired temperature.

V. SCOPE

The Installation and Operational Qualification protocol is intended to certify with documented evidence that the Drying Oven is installed properly and functions as desired by Giles.

VI. ROLES AND RESPONSIBILITIES



1. QA Laboratory

- ❖ Write and issue the protocol
- ❖ Investigate protocol deviation reports
- ❖ Execute the IQ and OQ.
- ❖ Review raw data and originate interim notification to Quality Assurance
- ❖ Write and route the final report

2. Quality Assurance

- ❖ Review and approve the protocol.

Controlled Document

	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 4 of 10	

- ❖ Review and approve raw data and notifications.
- ❖ Review, approve, and store the final report.

VII. TEST PROGRAM

A. INSTALLATION QUALIFICATION

Objective

The objective of the installation verification is to document the proper installation of the Drying Oven

Equipment/Materials

Quincy Labs 20GC Drying Oven

Procedure



Perform each item listed below for the Drying Oven.

- Location: Verify that the equipment is situated to allow sufficient room around the equipment for access doors and panels to be opened.
- Equipment: Document the Model and Serial or Asset Tag number of the Drying Oven.
- Utilities
 - Electrical Requirements: Verify the Drying Oven is receiving the specified voltage.

Acceptance Criteria

If the voltage is correct, each piece uniquely identified, and sufficient access for all doors and panels is available, the Drying Oven will be considered installed properly.

Controlled Document

	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 5 of 10	

B. OPERATION QUALIFICATION

Objective

The objective of Controls Verification is to document that the Drying Oven operates as needed by Giles. The controls will be operated to test the ability of the Drying Oven to be operated as needed.

Equipment/Materials

Quincy Labs 20GC Drying Oven

Procedure

Turn the Drying Oven on by pressing the power button.

Verify that the Drying Oven reaches desired temperature.

Acceptance Criteria

If the Drying Oven powers up and reaches the desired temperature then the Drying Oven is considered to be operationally qualified.

VIII. CALIBRATION

Verify that all instrumentation that requires calibration is calibrated.

Controlled Document

**GILES CHEMICAL ~ PREMIER MAGNESIA****Validation Protocol**

Title: Drying Oven IQ/OQ Protocol

Number: L13-VAL-100-054

Owner: Bryan Elchert

Revision: 0

Effective Date: July 17, 2013

Page: 6 of 10

PREMIER
MAGNESIA, LLC**Drying Oven: INSTALLATION QUALIFICATION****A. Installation Qualification****01. Location****a. Drying Oven:**

LOCATION			
Distance Criterion	Is the current area sufficient to open the access without obstructions (Yes/No)	Verified By	Date
Allow sufficient room around the oven for access doors and panels to be opened			

02. Equipment Identification

EQUIPMENT IDENTIFICATION			
Equipment	Serial or Tag Identifier	Verified By	Date
Drying Oven			
Comments:			

03. Utilities**a. Verify that the system is receiving its specified utility requirements.**

ELECTRICAL			
Specified	Actual	Verified By	Date
110 – 120 V Drying Oven			
Comments:			

Reviewed By: _____

Date: _____

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.

**GILES CHEMICAL ~ PREMIER MAGNESIA****Validation Protocol**

Title: Drying Oven IQ/OQ Protocol

Number: L13-VAL-100-054

Owner: Bryan Elchert

Revision: 0

Effective Date: July 17, 2013

Page: 7 of 10

**Drying Oven: OPERATIONAL QUALIFICATION****B. Operation Qualification****01. Controls Verification** -- to document that the Drying Oven controls work properly



CONTROLS/INDICATORS VERIFICATION				
Description	Function	Did Item function properly (Yes/No)	Verified By	Date
DRYING OVEN				
Power Button	With line power to the instrument, does pushing the Power Button cause the Drying Oven to heat up?			
Temperature Setting	With the Drying Oven powered on does it reach the desired temperature?			

Reviewed By: _____

Date: _____

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.

	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 8 of 10	



ATTACHMENT I - PROTOCOL DEVIATION REPORT LOG

Log each Protocol Deviation Report in the table below. Attach the PDRs to this Attachment.

PDR #	DESCRIPTION	DATE INITIATED	DATE RESOLVED
Comments:			

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.

	GILES CHEMICAL ~ PREMIER MAGNESIA		
	Validation Protocol		
	Title: Drying Oven IQ/OQ Protocol	Number: L13-VAL-100-054	
	Owner: Bryan Elchert	Revision: 0	
	Effective Date: July 17, 2013	Page: 9 of 10	

IX. PROTOCOL DEVIATION REPORT (PDR)

General Information

System Name: _____ Protocol Number: _____

Deviation Report Number: _____ Protocol Step & Page No.: _____

Instructions

Instructions The validation specialist assigns a sequential report number for each deviation with a specific protocol. For example, 001, 002, etc. can be easily referenced in a report.

1. Reference the relevant protocol number, step and page number of the noted deviation above.
2. Complete the below listed sections. If necessary, use additional pages and attach any supporting info.
3. Include the original PDR(s) with the protocol as an attachment. Summarize the impact of the deviation in the Validation Report.

=====

Description of Deviation:

Investigation Evaluation and Results:

Corrective Action and Resolution:

=====

Overall Investigation Review:

Prepared By: _____ Date: _____

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.



GILES CHEMICAL ~ PREMIER MAGNESIA

Validation Protocol

Title: Drying Oven IQ/OQ Protocol

Number: L13-VAL-100-054

Owner: Bryan Elchert

Revision: 0

Effective Date: July 17, 2013

Page: 10 of 10



ATTACHMENT III - SIGNATURE IDENTIFICATION LOG SHEET

Identify in the table below any personnel involved in the execution of this protocol.

Name	Affiliation	Signature	Initial	Date

Controlled Document

Only those quality documents viewed through the Giles Chemical electronic Documentation System are officially controlled. All other copies, whether viewed through another computer program or a printed version, are not controlled and, therefore, the Quality Unit at Giles assumes no responsibility for accuracy of the document.